Mesi ABI Doppler - The World's Fastest ABI Screening Device

Diagnose Peripheral Arterial Disease in **1** minute with MESI

Did you know you could make \$54.20 per test based on the Medicare item number 11610? Your practice could make \$10,840 per year based on 100 qualified patients being tested twice annually.

Ankle-Brachial Index (ABI) is a very simple comparison of blood pressure in the arm and legs. It is non-invasive and painless. With MESI ABI Doppler, the procedure becomes a reliable, objective screening device in primary healthcare.

Healthy arteries.

PAD is a circulatory problem where narrowed arteries reduce blood flow to your limbs.

MESI

Arteries are narrowed by build-up of plaque in the walls. The blood flow is partially restricted. The patient does not feel leg pain or other symptoms.

Arteries are clogged. The flow of oxygen rich blood is heavily restricted, leading to possible heart attack, stroke, gangrene, amputation and ulcerations.

PAD is likely to be a sign of a more widespread accumulation of fatty deposits in patient arteries (atherosclerosis). This condition may be reducing blood flow to their heart and brain, not only to their legs. The Mesi ABI Doppler can give you an early diagnosis, helping you determine the best treatment for your patient.

Simultaneous

measurement

Early diagnosis of PAD in primary healthcare is **cruicial**.



1 minute

measurements

No human

error



Healthcare staff

friendly



printout

Call us today to book in a trial and receive \$100 off the RRP!



Vital Medical Phone 1300 557 651 Fax 1300 557 631 sales@vitalmed.com.au www.vitalmedicalsupplies.com.au







MESI ABPI MD Screening Device

1 minute, non-invasive ABI screening for Peripheral Arterial Disease

Compared to the handheld Doppler probe, MESI ABPI MD is the world's fastest ABI screening device that performs an automated measurement. Empowered by a Plethysmography technique, it enables the device to provide accurate and objective results, based on which, physicians can diagnose Peripheral Arterial Disease with great confidence.

MESI ABPI MD KEY ADVANTAGES

Advanced error detection system

Smart software prevents false results even in the case of critical ischemia or medial calcinosis, giving physicians all the confidence they need.

Unique algorithm for ankle blood pressure calculation

It is not possible to measure blood pressure in ankles with a brachial blood pressure device. Therefore, the MESI ABPI MD algorithm is different, developed with human ankle anatomy in mind.

3Cuffs[™] - based technology


Plethismograpy sensors detect the smallest changes in volume. Ease-of-use excludes the possibility of human error.

Simultaneous measurement

Because blood pressure is constantly changing, simultaneous measurement is

crucial to avoid error; MESI ABPI controls the inflation and deflation of the cuffs at the exact same time to avoid blood pressure drift.



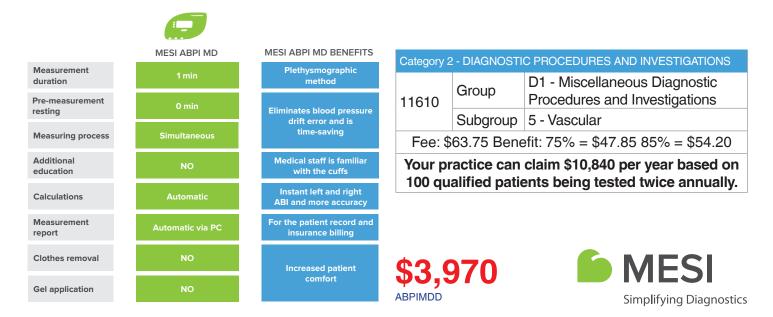

10% GST will apply to all items unless they are marked GST Free

Phone 1300 557 651 Fax 1300 557 631 sales@vitalmed.com.au www.vitalmedicalsupplies.com.au



Ability Dopplex

Ankle Brachial Index Measurements in approximately 5 minutes.



Dopplex[®] ATP Ankle & Toe Pressure Kit

The Ankle and Toe Pressure (ATP) Kit aids the assessment of arterial disease and neuropathy.

ATP Kit Includes:

- NEW DMX Digital Doppler
- Recharging kit
- EZ8-XS Widebeam 8MHz Doppler probe
- PPG Probe and adaptor
- Toe cuff and inflator
- Arm/Ankle cuffs (latex free)
- Trigger operated sphyg
- Neuropen and box of tips

3Year Warranty

... performance for life

HUNTLEIGH

Educational CD

\$3,450

• Gel

HUATPKIT

The Dopplex Bi-Directional Digital Doppler

- Integrated USB charging and connectivity
- · High sensitivity interchangeable probes
- Digital sound Processing
- High resolution colour display
- Bi-directional waveform generated from the digital Doppler spectrum with auto scaling feature
- Data storage on micro SD card







Hi-Bebe BT 200 Foetal Doppler

- LCD/OLED display
- High quality sound
- High sensitivity probe (2MHz)
- Hand-held style
- Low battery indicator
- Waterproof probe



Digital Handheld Dopplers – SR2 and SR3

Rate and Trace Display Dopplers with fixed 2MHz or 3MHz waterproof probe models.

These high sensitivity handheld Digital Dopplers are intended for general purpose screening of the fetal heart at all stages of pregnancy from as early as 8-10 weeks, right through to labour.

The attached probes are waterproof (IPX7) for high sensitivity FHR detection in both normal use and in waterbirth monitoring.

Choose between 2MHz (SR2) or 3MHz (SR3) fixed probe options. The 2MHz model is a general purpose probe for use at all stages of pregnancy. The 3MHz probe is the probe of choice for many clinicians for use in early gestation.

Features for both SR2 and SR3

- High Resolution Colour Display
- FHR BIG Numbers and Trace Display
- Attached waterproof probes
- Gel Suppression
- Record, Playback and Download of FHR & Sound
- Integrated Loudspeaker
- Sonicaid Software Compatible



SONICAID One HUNTLEIGH Rate display Doppler with fixed 2MHz probe

This compact foetal Doppler has been designed to fit easily into the clinician's pocket. Its wide beam coverage enables quick location of the foetal heart from early gestation to full term.

The SonicaidOne comes complete with soft pouch, gel, batteries, user guide, pregnancy calculator and pen.

Year

Warrant

www.huntleigh-diagnostics.com

- · Innovative and stylish
- Integrated probe and cable storage
- · Low running costs
- Pocket ergonomics

5625

HUSON1

· Large foetal heart rate display

Audio Doppler with fixed 2MHz waterproof probe



The traditional audio Doppler, setting the standard for routine intermittent auscultation at all stages of pregnancy.

Its high sensitivity probe and enhanced audio system provide optimum foetal heart detection, clear sounds and ease of use, from as early as 8-10 weeks gestation.

- Attached high sensitivity 2MHz probe
- Clear audio sound from integral loudspeaker
- Audio output for headphones

SONICAID D920

- Waterproof probe; ideal for waterbirth management
- · Carry case with gel, battery and user manual
- 5 minute auto-shut-off for maximising battery life
- 5 year warranty








HUNTLEIGH

HUNTLEIGH

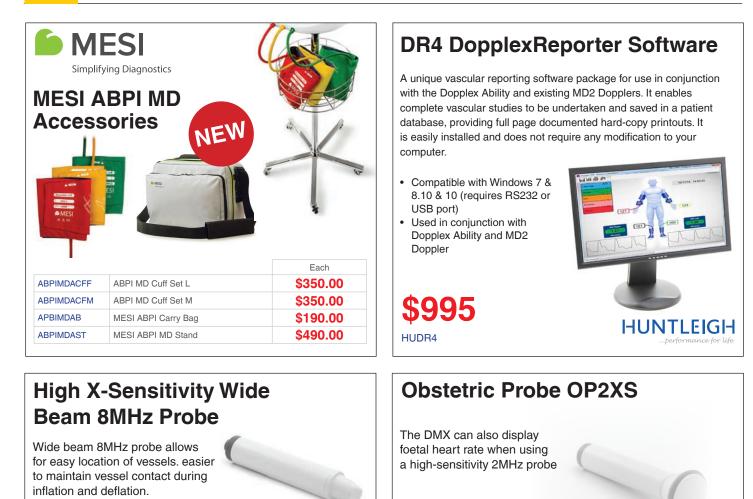

\$495 HUEZ8XS



Dopplex Ability Accessories

HUDASTAND	Stand and Trolley	\$695.00	HUDALAB	Ability Doppler Label Paper 5's	\$79.00
HUDALACSET	Ability Large Adult Cuff Set	\$575.00	HUDAPA	Ability Doppler Standard Paper 5's	\$55.00
HUDAFP	Fixing Plate	\$56.00	HUDASL	Ability Doppler Sleeves 100's	\$56.00

HLIOP2X9



Opening an account with us is easy.

Call **1300 557 651** and speak to one of our Customer Service team members. Alternatively you can open an account online at **www.vitalmedicalsupplies.com.au/register**