


Do you want a high quality, QCPP and Cold Chain compliant vaccine fridge at an affordable price?

Euro Chill are backed by 30 years experience in the industry so you can be confident that you will be buying the right Vaccine Fridge to fit your requirements and budget. Also known as Medical, Pharmacy or Immunisation Fridges, the Vacc-Safe® Premium line of Vaccine Fridges offer the following benefits:

- The best peace of mind warranty in the industry with 36 months parts and labour on the complete unit and a 5 year parts only warranty on the compressor
- The highest temperature and humidity rating of any vaccine fridge on the Australian market making it the ideal choice for the Australian climate
- Exceptional temperature stability with a “pipe & fin” cooling system instead of an older style “cooling plate” system found in most other brands of vaccine fridges
- Excellent security with a tight fitting secure lock
- QCPP and Cold Chain compliant
- Available sizes: 250 litre, 400 litre, 600 litre and 1000 litre

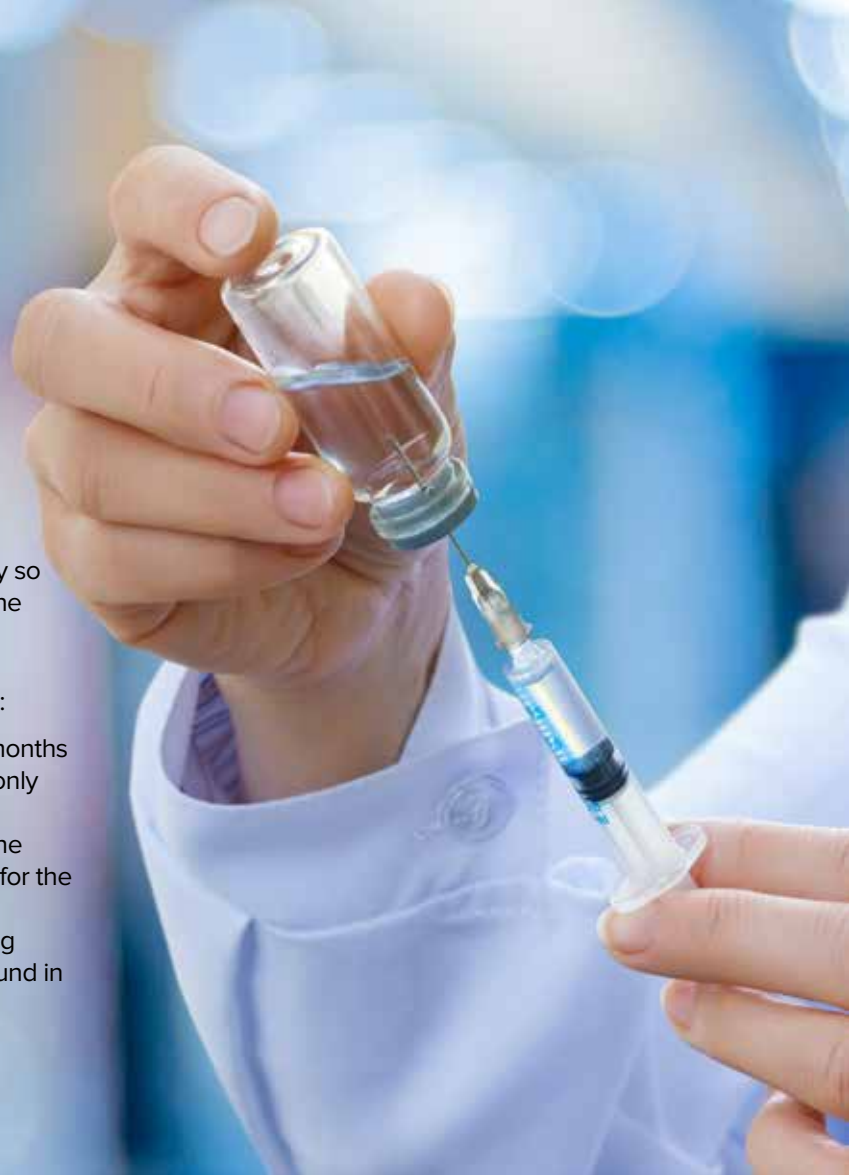

HIGH QUALITY STEEL (POWDER COATED) INNER LINER

Instead of the usual plastic liner, the Vacc-Safe® PREMIUM vaccine refrigerator comes with a powder coated steel liner. The steel liner is more durable as it is a better conductor of heat than plastic and also allows the vaccine containers to touch the liner without a loss of cooling efficiency.



BETTER LIGHTING

The Vacc-Safe® lighting has a shield on the lighting that directs the lighting inwards instead of outwards which gives exceptional illumination to every shelf without blinding the user.



IMPROVED SHELVING

The shelving is set up for small increments to allow for the best use of cabinet space. To stop vaccines from falling off, the shelves have an up stand at the front and rear. Optional shelf dividers available.



LARGE, EASY TO USE DIGITAL CONTROL

The Digital Controller is the largest and easiest to use on the market. New touch type buttons with a large temperature read out that is easy to see and use without straining your eyesight. It is flush mounted to the front of the cabinet so can be seen at any angle as opposed to other brands that recess their controller.



WELL-DESIGNED DOOR HANDLE

The door handle is ergonomically located at a convenient height and is fully moulded and concealed for functionality. Being concealed eliminates being caught by clothing or a lab coat and cannot be knocked or loosen over time.



EXTERNAL ALARM

Comes standard with a voltage free pair of contacts for external alarm. The contacts can be changed from NO to NC by pushing a button so no re-wiring required.



PORTABLE WITH FOUR CASTORS STANDARD

The Vacc-Safe® vaccine fridges all come with four (4) castors instead of the usual two (2) to allow for easy portable relocation and cleaning behind the unit.



OPTIONAL CLOUD-BASED WIRELESS LOGGING SYSTEM

This optional system gives you: Smart phone App, remote temperature monitoring with push notifications and/or Email Alert in case of Power Failure or Temperature Alarm Mode



Vacc-Safe® 250 Premium

- Gross volume 225L
- Dimensions W560 x D640 x H1451mm
- Operating range 2°C to 8°C
- 5 full size adjustable shelves
- High Ambient rating 40%/75% Humidity
- MEPS power consumption 1.428kWhr/24Hr
- Self-closing lockable doors
- Records Max/Min temperature history
- Audible and visual High/Low temp alarm
- Door open alarm
- Fan assisted cooling
- USB data logger and software included
- Supplied with NO/NC contacts for external alarm connection



NEW

Eur
Chill

\$3,300

VACCSAFE2

Freight charges will apply

Vacc-Safe® 400 Premium

- Gross Volume 400L
- Dimensions W560 x D640 x H1930mm
- Operating Range 2°C to 8°C
- 7 full size shelves
- High Ambient Rating 40°C / 75% Humidity
- MEPS power consumption 3.92kWhr/24Hr
- Self-closing lockable doors
- Records Max/Min temperature history
- Audible and visual High/Lowtemp alarm
- Door open alarm
- Fan assisted cooling
- USB data logger and software included
- Supplied with NO/NC contacts for external alarm connection



NEW

Eur
Chill

\$3,695

VACCSAFE4

Freight charges will apply

Vacc-Safe® 600 Premium

- Gross volume 600L
- Dimensions W700 x D710 x H1950mm
- Operating range 2°C to 8°C
- 6 full size shelves
- High ambient rating 40°C / 75% humidity
- MEPS power consumption 3.7kWhr/24Hr
- Self-closing lockable doors
- Records Max/Min temperature history
- Audible and visual High/Low temp alarm
- Door open alarm
- Fan assisted cooling
- USB data logger and software included
- Supplied with NO/NC contacts for external alarm connection



NEW

Eur
Chill

\$5,695

VACCSAFE6

Freight charges will apply

Vacc-Safe® 1000 Premium

- Gross volume 1000L
- Dimensions W1050 x D710 x H1950mm
- Operating range 2°C to 8°C
- 6 full size shelves per door
- High ambient rating 40%/75% humidity
- MEPS power consumption 7.4kWhr/24Hr
- Self-closing lockable doors
- Records Max/Min temperature history
- Audible and visual High/Low temp alarm
- Door open alarm
- Fan assisted cooling
- USB data logger and software included
- Supplied with NO/NC contacts for external alarm connection



NEW

Eur
Chill

\$7,500

VACCSAFE10

Freight charges will apply

BROMIC[®]

MediFridge™



Digital thermostat with min/max temperature recording



LED lighting

Free Data Logger

Fully configurable, multi-use temperature recorder that can create temperature reports



Visual & Audio Alarm



Self-closing door with lock



5 year warranty

Available in
145L, 220L &
374L
Capacity





Bromic Medifridge™ 374L + Data Logger

Bromics large capacity 374L MediFridge provides a solution for practices who require a certified medically graded refrigerator.

Features:

- Preconfigured temperature control system ensuring consistent temperatures are held on all shelves between 2°C & 8°C, which is a QCPP approved specification
- Open rung style shelving units allows for greater circulation of cold air
- Fan forced air increases the circulation around the shelves
- Audible alarms when temperature operates out of predetermined range
- Adjustable shelving spaces
- External LCD shows the minimum and maximum readings
- Self closing lockable door
- Right hand hinged glass door for easy viewing
- Internal light
- Additional Information:
- Width: 595mm
- Depth: 598mm
- Height: 1840mm



\$3,295

BR374GD

Freight charges will apply

Bromic MediFridge™ Display Chiller 145L

- Integrated temperature Alarm System ensures sensitive stock is maintained and managed according to requirements
- Fan Forced air circulation for even temperature throughout the unit
- Increased air circulation around the unit thanks to unique custom wire 'open rung' shelves



\$1,865

BR140GD

Freight charges will apply



Call us on 1300 557 651

Our Customer Service team are available to assist you with orders, enquiries and quotes from 8.30am to 5pm, Monday to Friday (excluding Public Holidays).



Medi-Log II Logger Thermometer

- Medi Log II - Automatic setting data logger for Vaccines and Food.
- Meets all the latest Strive for Five guidelines and WHO requirements
- Fully self contained, sets its own parameters automatically for Vaccines
- One button to use for all functions
- No extra cradle needed.
- Free software
- Includes Glycol Vial, 2 Meter External probe and USB Cable



\$187.95

HLMEDILOG2

Large Display timer



\$10.45

HLBIGTIMER

Medi-Temp Fridge/Freezer Digital Thermometer

- 2 temperature sensors: 1 inbuilt inside the unit and another sensor on a 3 meter wire.
- Inbuilt alarms that can tell you that the temp limits are exceeded.
- Inbuilt Memory for temperature history readings.
- Inbuilt sensor range measures from: -10 up to +50C and,
- Outer sensor (on the end of the cable) measures temperatures from -50 C up to +70C.



\$31.90

HLT009

Wireless Temperature Data Logger



\$385

HLENLESS30

Make sure you meet your Strive for 5 requirements

HLP Controls
Innovating solutions in testing technology

