



Choose the right stethoscope to help you perform at your best.

3M™ Littmann® Stethoscopes feature adjustable headsets with soft-sealing eartips. Select Littmann mechanical stethoscopes are designed with tunable technology to give you the ability to hear different frequency sounds by simply adjusting the pressure applied to the chestpiece. Get the acoustic performance, high quality and comfort you need with a Littmann stethoscope.

Acoustic Performance Levels



Classic II Pediatric and Infant Stethoscopes	Master Classic II™ Stethoscope	Classic III™ Monitoring Stethoscope	Cardiology IV™ Diagnostic Stethoscope	Master Cardiology™ Stethoscope	Electronic Stethoscope Model 3200
Hard-working, dependable, clinical tools for physical assessment and monitoring of your patients, one long shift after another.			Provides outstanding acoustic sensitivity to help you perform detailed diagnostic auscultation. Useful in challenging and critical care environments.		High-performance auscultation. Hear more of what you need to hear with amplification of 24x, ambient noise reduction and optional connectivity.

Approximate Length	71cm	68cm	69cm	69cm	68cm	68cm
Approximate Weight	105g, 95g	160g	150g	177g	185g	185g
Warranty	3 years	3 years	5 years	7 years	7 years	2 years
Diaphragm	Non-Tunable (Floating)	Tunable	Tunable (2)	Tunable (2)	Tunable	Electronic
Approximate Chestpiece Diameter	Double-sided Paediatric 3.4cm, 2.7cm Infant 2.8cm, 1.7cm	Single-sided 4.3cm	Double-sided Adult side 4.3cm Pediatric side 3.4cm	Double-sided Adult side 4.3cm Pediatric side 3.4cm	Single-sided 4.8cm	Single-sided 5.1cm
Features	Specially designed for the smallest patients	Designed for easier grip and maneuvering	Convertible open bell	Convertible open bell Dual-lumen tubing	Dual-lumen tubing Best acoustic response in Littmann mechanical line	Eliminates 85% (on average) of ambient noise ¹ Bluetooth® wireless technology Recording and playback



The power to learn, in the palm of your hand.

3M™ Littmann® Learning Institute App – Free with purchase of a 3M™ Littmann® Stethoscope!

1. 3M data on file. 3M Australia Pty. Limited.

Need more help? Visit Littmann.com.au/FindMyStethoscope



3M™ Littmann® Cardiology IV™ Stethoscope

The 3M™ Littmann® Cardiology IV™ Stethoscope is the next generation of one of the most recognized stethoscopes in the profession, offering improved acoustic performance with better audibility of high-frequency sounds (adult side, firm pressure).

Its distinctive chalice-shaped chestpiece features tunable diaphragms on both sides to deliver outstanding acoustic performance, diagnostic versatility, and convenience for treating adult and paediatric patients. The paediatric side of the chestpiece converts to a traditional open bell by removing the tunable diaphragm and replacing it with the included non-chill bell sleeve.

Features:

- Outstanding acoustics with better audibility of high-frequency sounds (adult side, firm pressure).
- Tunable diaphragm: Hear high or low frequency sounds by slightly adjusting pressure on the chestpiece.
- Dual-lumen tubing: Two sound paths in one tube eliminates the rubbing noise of traditional double tubes.
- Next-generation tubing lasts longer and resists stains. Not made with natural rubber latex or phthalate plasticizers.
- Headset is easily adjusted for individual fit and comfort. Angled eartubes align with ear canals.
- 7 year warranty



		Each
3M6152	Black	\$311
3M6154	Navy	
3M6155	Hunter Green	
3M6156	Plum	
3M6158	Raspberry	



3M™ Littmann® Classic III™ Stethoscope

The Littmann Classic III stethoscope is the latest version of the stethoscope that helps millions of medical professionals achieve their best. It brings new design, materials, and technology to the Littmann stethoscope series trusted by clinicians around the world.

The Classic III stethoscope offers you:

- High acoustic sensitivity
- Tunable diaphragms
- Traditional & Paediatric Bell options
- Latex Free
- Next-generation tubing
- Quality & reliability
- Available in multiple tubing colours and chestpiece finishes

		Each
3M5620	Black	\$152.43
3M5622	Navy	
3M5627	Burgundy	
3M5832	Lavendar	
3M5835	Turquoise	
3M5839	Lemon	



3M™ Littmann® Master Cardiology™ Stethoscope

- Two-in-one tube design (dual lumen) eliminates noise from two tubes rubbing together
- Pressure-based sound frequency adjustment with tunable diaphragm
- Anatomically designed headset
- 7 year warranty

		Each
3M2160	Black	\$311
3M2163	Burgundy	
3M2164	Navy Blue	





3M™ Littmann® Electronic Stethoscope



Model 3200



State-of-the-Art Sound Sensor

Providing a life-like sound experience equivalent to a high-end cardiology stethoscope



Frictional Noise Reduction Technology

Designed to reduce distraction from handling noises



Up to 24X Sound Amplification

Helpful for faint heart, lung or body sound, obese patients or when a patient's clothing restricts your listening experience



Proprietary Ambient Noise Reduction Technology

Reducing on average, 85% of interfering background noise

- The 3200 delivers exceptional sound quality and a simple-to-operate user interface. Clinical evidence shows that it's easier to detect difficult-to-hear heart sounds, like S3 gallops, aortic regurgitation murmurs and abnormal lung sounds. With these sound and usability improvements, you might never go back to your old stethoscope.

Features of the 3200 include:

- On-board recording, now with 12 tracks and playback capabilities Helpful when seeking a peer opinion or when you are writing up notes on your patient's condition, after your patient has left
- Bluetooth® technology lets you transmit sounds real-time to your PC, which can then be further analysed or attached to medical records



\$684

3M3200BK




3M™ Littmann® Classic II™ Paediatric Stethoscope

The 3M Classic II Paediatric stethoscope is a model that has a two sided chestpiece with the traditional bell and diaphragm and is designed and sized to provide Littmann quality acoustics for the pediatric population. This stethoscope can be used for physical assessment and diagnosis.

Features

- Two sided chestpiece with the traditional bell and diaphragm
- The large diaphragm has a 1.31 inch or 3.3cm diameter. The bell or small diaphragm has a diameter of 1 inch or 2.5cm.
- Snap tight soft-sealing eartips create an excellent acoustic seal to reduce ambient noise and provide a comfortable fit.
- Nonchill rim and diaphragm for the comfort of your patients.
- Weight: 3.70oz

		Each
3M2113	Black	
3M2119	Caribbean Blue	\$152.39
3M2122	Raspberry	

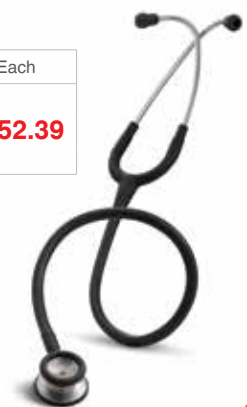

HEINE GAMMA C3 Cardio Stethoscope

Exceptional acoustic performance, particularly in the cardiac frequency range.

- Combined double chest-piece with two different frequency ranges.
- Instant selection of a double membrane on one side or a conical acoustic head on the other.
- High-performance double membrane, 45 mm diameter gives ideal sound transmission and perfect skin contact for exceptional sound reproduction.
- Two-in-one tube system for enhanced performance. Very supple, yet kink-resistant with exceptional acoustic reproduction.
- Includes ear olives (large / soft) and one pair of replacements (small / hard) for sound insulation.



\$359

HEM0009944





HEINE GAMMA® 3.1 Pulse Stethoscope

Stethoscope with solid flat chest piece. The chest piece can be easily inserted, under the sphygmomanometer cuff.

- High-quality membrane 44 mm diameter
- Single-tube construction with excellent sound transmission
- Includes ear olives (large / soft) for sound insulation



\$98

HEM0009941



HEINE GAMMA® 3.2 Acoustic Stethoscope

Combined double chest piece with two different frequency ranges to choose from:

- A high-quality membrane on one side and a conical acoustic head on the other side.
- Single-tube construction with excellent sound transmission.
- High-quality membrane, 44 mm diameter for excellent acoustic performance and sound transmission.
- Cone diameter 30mm.



\$127

HEM0009942



HEINE GAMMA® 3.3 Paediatric Stethoscope

- Combined double chest piece with two different frequency ranges to choose from:
- A high-quality membrane on one side and a conical acoustic head on the other side.
- High-quality membrane, 31 mm diameter for excellent acoustic performance and sound transmission.
- Cone diameter 30mm.



\$127

HEM0009943



Sprague Rappaport Stethoscope

- High, medium and low frequency detection
- Complete with pouch, 3 bells, 2 diaphragms and 2 pairs of ear tips
- Latex-free tubing
- More colours available



SA0418	Black
SA0418BY	Burgundy
SA0418B	Navy Blue
SA0418GY	Grey

\$39



Classic Stethoscope

- Stainless steel chest piece
- Double inner spring binaural
- Spare ear tips

ADULT
Black

\$59

SA0417



Cardiology Stethoscope

- Extra deep bell complete with non-chill rings on both sides
- Super sensitive diaphragm for low frequency sounds

SA0421	Black
SA0421BY	Burgundy
SA0421G	Green

\$99





Harvey Elite Stethoscope

The Harvey Elite offers updated styling with the superior comfort and auscultation capabilities cardiologists have come to rely on. Its stainless steel head transmits sounds crisply and clearly.

- Full-range cardiology/pulmonology stethoscope
- Double head (bell and diaphragm)
- Paediatric chestpiece included
- Stainless-steel chestpiece with non-chill rim
- Rotatable binaurals and interchangeable comfort sealing eartips
- Dual-lumen tubing creates individual, uninterrupted sound channels to each ear
- Flat or corrugated diaphragms
- Includes 3 extra pairs of ear tips
- Latex free
- Heart sound CD included
- 71cm length



10 Year Warranty

	Each	3+
WL5079125S	\$226	\$209



Training Stethoscope with twin binaurals



\$35

SA0565



Dual Head Stethoscope

		Each
SA0416	Black	\$22
SA0416B	Blue	
SA0416R	Red	

LATEX FREE

Dual Head



Arden Double Head Stethoscope

- General-use stethoscope
- Lightweight chest pieces with nonchill rims
- Latex-free, flexible, single-channel tubing
- 1 year warranty
- 63.5cm length



WLARDEN-BK	Black	\$22.05
WLARDEN-BL	Blue	



Pinard Foetal Stethoscope

SAS420	Plastic	\$7.95
SAS426	Aluminium	\$9.95



Nurse Stethoscope

		Each
SA0415	Black	\$12.95
SA0415B	Blue	
SA0415Y	Yellow	
SA0415R	Red	
SA0415G	Green	



Lightweight Stethoscope

- General-use stethoscope
- Lightweight chestpieces with nonchill rims
- Adjustable chrome-plated binaurals
- Latex-free, flexible, single-channel tubing
- 81cm length

WL5079-73	Velvet Black	\$37.80
WL5079-76	Misty Blue	
WL5079-74	Poppy Red	



2 Year Warranty



Professional Stethoscope

- General-use stethoscope
- Stainless steel double-head (bell and flat diaphragm) chestpiece with nonchill rim
- Rotatable binaurals and interchangeable comfort sealing ear tips
- Latex free
- 5 year warranty
- 71cm length



WL5079-135	Black	\$103.95
WL5079-137	Navy	\$103.95

