

HEINE NC2 Dermatoscope

All-round dermatoscope, universal application, efficient and precise in all situations

Experience The Difference



LED HO

LED NOW IN HEINE QUALITY.

- Fast change from contact to non-contact examination
- Perfect combination of precise optics and bright illumination
- Available with digital documentation





NC2 Dermatoscope

The most versatile Dermatoscope available

All-round dermatoscope, universal application, efficient and precise in all situations: The contact plate achieves a clear, reflection-free and highly detailed image with 10x magnification, even without immersion fluid. Without the magnetically-held contact plate, the NC2 is used for non-contact examinations in 6x magnification for an efficient examination with a large field of view of the skin. The iPhone* accessory for the Apple iPhone from Version 5* and the NC2 app enable extensive photo documentation, in non-contact mode too.

- Fast change from contact to non-contact examination
- Perfect combination of precise optics and bright illumination
- Available with digital documentation

HED8417812	NC2 Deramatoscope with C/Plate	\$1,435.00
HED8427812	NC2 Deramatoscope without C/Plate	\$1,330.00



LEDHQ
LED NOW IN HEINE QUALITY.



iC1 Dermatoscope and Documentation System

PERFORMANCE:

- > View or show your patients captured images directly on screen.
- > HEINE precision optics. High-quality aspherical optical system with up to 40 x magnification and up to 12 mega pixel image quality.
- > LEDHQ Illumination. 4 high-performance-LEDs for homogeneous, bright illumination and excellent image: Colour temperature typ. 5,000 K, typ. CRI > 80.

FUNCTION:

- > Free Onboard App: HEINE iC 1
 - Easy location of images on Body Map
 - Input of patient information
 - E-mail function for sharing images
 - Adding of comments and notes to images
 - In-app magnification
- > Changeable and cleanable contact plate made of high quality glass. With a scale for exact image comparison.
- > Future compatibility – with exchangeable mobile phone cases (currently available for iPhones 5 / 5s / 6 / 6s / 7 / SE / 8).
- > Toggle Function to conveniently alternate between polarised and non-polarised illumination.
- > Latest Li-Polymer battery technology with USB charging function and charge status indicator.



LEDHQ
LED NOW IN HEINE QUALITY.



HEK2702830	iC1 Dermatology Set (6/6S)	\$1,435.00
HEK2712830	iC1 Dermatology Set (5/5S)	
HEK2732830	iC1 Digital Dermatology Set (7/8)	

Mini 3000 LED Dermatoscope

Features:


- Solid metal construction of the instrument head. Durable and long lasting.
- High quality, achromatic optical system with 10x magnification.
- Choice of two different standard contact plates.
- Each mini 3000 Dermatoscope includes a comprehensive HEINE dermoscopy compendium with explanations and photos.
- Twice as bright as a conventional XHL Xenon-Halogen instrument.
- True and natural colour rendering: Red is red, blue is blue.

Contents

- Contact plate
- Battery handle (for AA batteries)
- Dermatoscopy Oil
- Zipper case

\$715

HED8778021

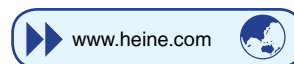

LEDHQ
LED NOW IN HEINE QUALITY.





Delta 20® T Dermatoscope

Highest Quality Images and Diagnostic Flexibility

PERFORMANCE:

- **Precision HEINE Optics** – High quality, achromatic optical system provides unmatched sharpness and resolution. Focusing system provides for 10 x to 16 x magnification for recognition of even the smallest details and structures.
- **LED in HQ** – The new standard in LED illumination defining optimal light intensity, homogeneity and colour rendering for each individual examination situation – Colour temperature: 5.000 K, Colour Rendering Index (CRI) > 80 on a scale of 100. This ensures true tissue colour for the most accurate diagnosis.

FUNCTION:

- **Toggle between Polarised and non-Polarised illumination modes** – Conveniently alternate between Polarised and non-Polarised views with the press of a button – no exchange of contact plates necessary. Allows recognition of the “Blink Sign” when observing crystalline structures and presence of milia cysts, increasing diagnostic function.
- **Optional small contact plate available** – Special small contact plate for difficult to access lesions.
- **Digital documentation capabilities** – Connect to major brand SLR digital cameras with HEINE digital camera lens system. Specially designed lenses integrate the camera and DELTA 20 T optics together to provide the highest quality imaging possible.

DURABILITY:

- **HEINE offers a 5 year guarantee**
HEINE makes no compromises in manufacturing high quality medical instruments. Our commitment to vertical integration in manufacturing means we control all aspects of our instrument quality, from utilizing carefully selected, well-matched materials to a high level of hand assembly. This ensures that each HEINE instrument meets or exceeds all requirements in any medical environment.
- **LED Thermal Management** – LED illumination systems require thermal management to ensure long-term performance and product stability. Carefully selected and applied heat-diffusing materials in the DELTA 20 T guarantee consistent quality and performance throughout the LED's lifetime rating of up to 50.000 hours.



LED HQ

LED NOW IN HEINE QUALITY.

DELTA 20 T Dermatoscopy Set

Includes:

- HEINE DELTA 20 T Dermatoscope Head
- Contact plate with scale
- BETA battery handle
- Dermatoscopy oil 10ml
- Hard case

\$2,095

HEK6210118



DELTA 20 T Rechargeable Dermatoscopy Set

Includes:

- HEINE DELTA 20 T Dermatoscope Head
- Contact plate with scale
- BETA 4 NT rechargeable handle
- NT 4 table charger
- Dermatoscopy Oil 10 ml
- Hard case

\$2,795

HEK6224420



DELTA 20 T USB Rechargeable Set

Includes:

- HEINE DELTA 20 T Dermatoscope Head
- Contact plate with scale
- BETA 4 USB rechargeable handle
- USB Cord and plug In power puppy
- Dermatoscopy Oil 10 ml
- Hard case

\$2,445

HEK6228388



DELTA 20 T Dermatoscopy Head Only



HEK0834221

Contact Plate with Scale

\$1,835



3Gen



Polarized or non-polarised?

Cross-polarisation cancels out the light that is being reflected by the skin surface, thus allowing a deep, glare-free view into the skin—without immersion fluids.

Non- polarised light requires a glass contact plate & immersion fluids to see lesion structures.

Leading dermatologists and Skin Cancer Specialists increasingly recommend devices that allow for both cross-polarised and non-polarised dermoscopy, and studies are showing diagnostic advantages from having both modes available. In practice however, you'll most often stick to polarized light.

Camera compatibility

Having the ability to take a photo of a special case (prior to excision), with a camera, smartphone or iPad®, is powerful. DL100 or Carbon cannot be attached to those devices.

Magneti-Connect

DermLites are easily able to magnetically connect to any Camera, iPhones, Samsung Galaxy Phones and even iPads.

- Made in the USA
- 5 Year Warranty
- Local Australian Service Centre



DermLite DL100

White LED lighting technology combined with a 10x Hastings Triplet lens provides for ideal skin microscopy. The LED daylight-type light is linearly Polarised, and a second polarizing lens is aligned to cancel the reflected light from the skin. A high quality, three-element lens with 10x magnification, superb color correction, and reduced image distortion is used to image the skin.

\$389

EQGENDL100

dermlite®



DermLite DL3N

The DL3N features a 25 mm four-element lens, 28 high-powered LEDs and the all-new PigmentBoost™ illumination. For optimum ease of use, a quick tap on its power button instantly toggles between polarization modes while the optional PigmentBoost mode produces a warmer illumination reminiscent of conventional dermatoscopes.

\$1,545

EQDL3N

dermlite®



DermLite Foto II Pro

Featuring Polarised and non-Polarised illumination, a 6-element lens design for superior resolution, zoom, camera-controlled automatic polarization change, Foto II Pro is the go-to device for situations where only the best dermoscopic image quality will suffice.

\$2,750

EQDERMLI29

dermlite®





DermLite DL200 HR

DermLite DL200 HR is permanently cross-polarized, with two brightness levels, making it the brightest polarized DL200.

Key benefits:

- Compact, ergonomic design
- Mobile device-compatible
- >20% more life on a charge than DL2
- Four-level charge indicator
- Automatic shutoff to save battery life
- Contact & non-contact dermatology
- Faceplate with 10 mm markings
- Ultra-bright & polarized only

\$1,250

EQDL200HR



dermlite®

DermLite DL200 Hybrid

DL200 Hybrid combines polarized & non-polarized lighting. Visualize pigment and vascular structures with or without immersion fluids and see the benefits of being able to toggle between polarized and non-polarized illumination.

Key benefits:

- Compact, ergonomic design
- Mobile device-compatible
- >20% more life on a charge than DL2
- Four-level charge indicator
- Automatic shutoff to save battery life
- Contact & non-contact dermatology
- Faceplate with 10 mm markings
- Ultra-bright & polarized only

\$1,375

EQDL200HY



dermlite®

DermLite DL4

- Compact all-aluminum design
- Now 45% brighter (DL4W only)
- Instant on and auto shut off
- Large 30mm 10x lens system
- Charge from any USB Port
- Polarized and non-polarised
- PigmentBoost™ Plus (not in DL4W)
- Contact and non-contact dermatology
- Magnetically attach cameras, smartphones and iPad
- IceCap™ infection control system
- Four charge indicator levels
- Two brightness levels
- Magnetic faceplate
- Sturdy protective belt holster

Versatile DL4 vs Ultra-Bright DL4W

While model no. DL4W produces 45% more light from its exclusively white LEDs, the top-of-the-line model DL4 provides enhanced imaging of deeper pigmentation with two levels of colour spectrum control (PigmentBoost™ Plus).



EQDL4	DL4	\$1,959
EQDL4W	DL4W	\$1,795

dermlite®

DermLite Lumio

Designed to view the skin with magnification and clarity. A ring made up of 40 LEDs and polarization filters provide glarefree, cross-Polarised light by canceling the reflected light from the skin. A high quality, 75 mm 2x lens with superb colour correction and reduced image distortion.

EQLUM	Lumio	\$675
EQLUMUV	Lumio UV	

dermlite®





Heine iC1 Accessories



HEX0099303	E4-USB Power Supply + USB Cable	\$86
HEK0034201	Contact Plate Heine iC1 with scale	\$315
HEK0034250	Mobile phone case for iPhone 6 Series	\$86
HEK0034251	Mobile phone case for iPhone 5 Series	
HEK0034252	iPod Touch 6 Case	
HEK0034253	Mobile phone case for iPhone 7 Series	



Delta 20 T Accessories



HEK0034205	Contact plate small, 8mm dia. (for the examination of inaccessible lesions)	\$199
HEX0099008	BETA Belt Clip	\$52
HEX0009923	Adaptor Cord (1.3m)	\$84
HEK0034186	Nikon SLR Photoadaptor for Delta 20 Plus	\$765
HEK0034185	Canon SLR Photoadaptor Delta 20 Plus	\$765



Dermatoscopy Oil

Dermatology oil is applied to external skin layers, making them transparent.



Each
\$17.95
HEK0034005



Dermlite Accessories



EQDLCONK	Connect Kit iPhone 4	\$64
EQDLCONK5	Connect Kit iPhone 5	
DERMLI-24	Connect Kit iPhone 5C	
EQDLCKI6	Connect Kit iPhone 6	
EQDERMLI33	Connect Kit iPhone 6 Plus Adaptor	
DERMLI-66	Connect Kit iPhone 7	
DERMLI-67	Dermlite 7 Plus Phone Case	
DERMLI-72	Connect Kit iPhone X	\$80
DERMLI-18	Connect Kit Galaxy S3	
DERMLI-22	Connect Kit Galaxy S4	
DERMLI-31	Connect Kit Galaxy S5	
DERMLI-58	Connect Kit Galaxy S6	
DERMLI-64	Connect Kit Galaxy S7	
DERMLI-50	Connect Kit Galaxy S8	
DERMLI-55	Connect Kit Galaxy S8+	
DERMLI-20	Connect Kit iPad MINI Adaptor	
DERMLI-60	Connect Kit iPad MINI 4	
DERMLI-19	Connect Kit iPad 2/3/4	\$110
DERMLI-25	Connect Kit iPad AIR Adaptor	
DERMLI-38	Connect Kit iPad AIR 2 / PRO	
DERMLI-52	Connect Kit iPad AIR PRO 10.5"	
DERMLI-53	Connect Kit iPad AIR PRO 9.5.5"	
DERMLI-61	Connect Kit iPhone Touch	\$64
EQDERMLI59	Icecap Infection Control Caps for DL4 (Box/25)	\$22
DERMLI45	Universal Smart Phone Adaptor	\$120
CAMERA-02	Universal Camera Adaptor Kit	\$95

dermlite®

Download the FREE Vital Medical Supplies APP
VISIT PAGE 10 FOR MORE INFORMATION

NEW
FREE download available from:

Download on the App Store
GET IT ON Google Play