

Medisafe **PLUS**

WILL KEEP YOU SAFE ...!



MEDISAFE PLUS is a purpose-built fridge, designed to Australian standards and specifically for the storage of vaccine and medical products. Our Medisafe Plus range are feature packed, resulting in a vaccine storage solution that gives you peace of mind, knowing your vaccines are stored in a controlled environment. Medisafe Plus units are available in different sizes to suit your all your storage solutions.

“Strive for Five (5°C)”

Vaccines lose their efficacy or may be destroyed if they are frozen; get too hot or are exposed to direct sunlight or fluorescent light. Vaccines must be stored and transported within the recommended temperature range of +2°C to +8°C.

“Strive for 5°C” as this allows the greatest leeway in temperature fluctuation.

The Medisafe Plus system

Our Medisafe Plus refrigeration equipment is installed with two electronic controllers each ensuring temperature in the cabinet is controlled and monitored. Air temperature can change quickly in response to the opening of the door or activation of the refrigeration system. However, air temperature is NOT a measure of the core temperature of the vaccine, this can vary pending on the environment temperature. One controller uses these measurements to cycle the refrigeration system to maintain the temperature within the required parameters. A suitable delay is programmed into the system to prevent unnecessary false alarms due to opening of the door or other legitimate events. The other controller is a secondary visual monitor of the refrigeration air temperature, it also assists in visual reading of the Min and Max temperature of the unit

AQ Box

The AQ Box is a device that is installed to a Medisafe Plus refrigeration unit and enables a user to have remote information sent to both monitor the cabinet temperature and also send notifications of faults with in the cabinet. This system can monitor the unit and you will be informed via an email or text, if there is a power failure, or your product / cabinet temperature is fluctuating outside the set temperature parameters. It automatically logs data to the Cloud where you can access information anywhere you have internet access. With minimal monthly cost the AQ Box enables you to respond quickly to ensure product integrity.



Features of the AQ box include:

- Email and Text message alarm system to immediately notification staff of temperature breaches and power outages.
- 2 x RTD Class A temperature probes measuring both air temperature and simulated product temperature
- 8 hour Battery backup activated during power outages.
- Cloud based reporting accessibility by any computer, tablet or smartphone in any location.

* Requires a stable 3G mobile data signal to operate AQ box cloud monitoring



Compliant Pharmacy Refrigeration System

Features

- Superior Quality; No domestic fridge components used
- Electronic temperature control and monitoring system ensures consistent temperatures are maintained on all shelves
- Built-In cloud based temperature monitoring system (system contains its own independent wireless connection and does not rely on stores having their own WIFI connection)
- Provides continuous temperature monitoring.
- Email and SMS alerts notifying of temperature breach and power outage
- Min/Max Temperature Recorder logger for ease of daily temperature recordings Visual and audible auto resetting alarms and a power failure alarm
- VI311MS and VI371MS models include an internal evaporator guard and door fan switch
- All fridges come standard with an integral lock and key
- 24 hour technical support & service
- 24 months warranty
- Australian design *Manual Data Logger option available upon request *Subscription to cloud based reporting system required (min 3 year term)

NEW


Bigger is Better!

1500L

\$9,995

VI1365MS

Height	210cm
Width	150cm
Depth	117cm
Litres	1500

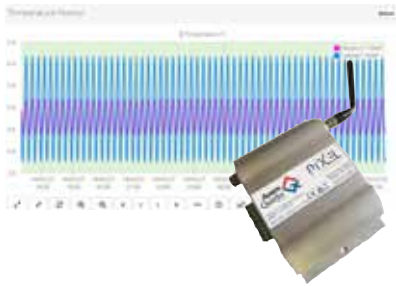

CLOUD BASED REPORTING INSTALLED
 *Requires a stable 3G mobile data signal to operate AQ box cloud monitoring

80L

\$2,595

VI80MS

Height	83cm
Width	50cm
Depth	56cm
Litres	80

306L

\$3,395

VI311MS

Height	163cm
Width	59.5cm
Depth	59.5cm
Litres	306

381L

\$3,695

VI371MS

Height	193cm
Width	59.5cm
Depth	59.5cm
Litres	381

NEW

700L

\$6,995

VI625MS

Height	210cm
Width	75cm
Depth	82cm
Litres	700

Auto Dialler (optional) \$490 VIDIALLER

Freight charges will apply

When the controller detects a critical condition the auto-dialler will be activated. This remote monitoring device will automatically call you to notify loss of power. It will dial up to five different pre-programmed numbers until the alarm is acknowledged. The auto-dialler requires a standard analogue phone to operate, and can be connected in parallel with other devices. It does not rely on the phone line for power but uses standard AA batteries.

10% GST will apply to all items unless they are marked GST Free

Phone 1300 557 651 Fax 1300 557 631 sales@vitalmed.com.au www.vitalmedicalonline.com.au





Strive for 5



Wi-Fi Temperature Sensor Data Logger with Glycol Vial
Model: B10-200

- Purpose built temperature monitor for vaccine and medical storage
- Connects to existing WiFi and logs all data to a cloud server.
- Alarms, visual, audible and email.
- Complies with all requirements for medical / vaccine storage.
- Data and graphs are available from anywhere with internet access, phone, tablet or computer.
- FDA21CFR, NIST, ISO17025
- Currently used in more than 1,000 hospitals.

B10 Wireless Temp Monitor

Purpose Built temperature monitor for Vaccine and Medical Storage.

- Connects to Existing WiFi.
- Logs all Data to Cloud.
- Alarms, Visual, Audible and Email.
- Complies with All requirements for Medical and Vaccine Storage.

\$675

HLP-B10



NEW

B11 Wireless Temp Monitor

Features

- Dual Point WiFi Logger.
- Logs all Data to Cloud. -
- 40 to +100C.
- With Display and Alarms.
- Mains and Battery.

\$850

HLP-B11



NEW



Medi-Temp Fridge/Freezer Digital Thermometer

Features

- 2 temperature sensors: 1 inbuilt inside the unit and another sensor on a 3 meter wire.
- Inbuilt alarms that can tell you that the temp limits are exceeded.
- Inbuilt Memory so that you can see what the temperature got up to and down to during the weekend or overnight.
- Inbuilt sensor range measures from: -10 up to +50C and,
- Outer sensor (on the end of the cable) measures temperatures from -50 C up to +70C.
- With magnetic and velcro mounting.

NEW

\$30

HLT009



Medilog Fridge/Freezer Thermometer

Features

- Pre set alarm set points for the vaccine requirements (alarm set points can be changed by the user).
- Shows Minimum and Maximum temps reached and a display of the time that the temperature has been outside the set area.
- Includes a 3 meter long small wire probe that can go through the fridge seals, a small Velcro probe holder, and a stick on probe holder. Temp range -20C to + 70C. (C/F Selectable).
- Battery powered - Immune to power outs.

NEW

\$75

HLMLOG



PDF Data Logger

Features

- Self Contained Pen Style data Logger.
- 24000 Readings.
- -50 - +70 C. Min/Max Readings.

NEW



AZ88161	With External Probe	\$130.00
AZ88160	Without External Probe	\$125.00

Temperature and Humidity Logger

Features

- Temperature and RH Data Logger
- Remote probe.
- AC240V Charger and Battery Back Up.
- SD Memory Card included.

NEW



\$195

HLP87799



New Facility or upgrading your practice?

We can help you make the right decisions and bundle a special equipment package for you at a great price!

Call 1300 557 651