

ORDERING INFORMATION

Leukosan® Strip White	Size (mm)	Strips/Pouch	Pouches/Box	Boxes/Shipper
72629-00	6 x 38	6	50	4
72629-01	3 x 75	5	50	4
72629-02	6 x 75	3	50	4
72629-03	6 x 100	10	50	4
72629-04	12 x 100	6	50	4
72629-05	25 x 100	4	25	4



Leukosan® Strip Tan	Size (mm)	Strips/Pouch	Pouches/Box	Boxes/Shipper
72629-07	6 x 38	6	50	4
72629-08	6 x 75	3	50	4
72629-09	12 x 100	6	50	4



Leukosan®

Wound Closure

New
Leukosan® Strip

The new generation of wound closure strips



BSN medical Australia & New Zealand

L6, 211 Wellington Road,
Mulgrave, Victoria, Australia 3170

www.bsnmedical.com.au
www.bsnmedical.co.nz

BSN0065

05/2016

Customer Service Australia
T 1300 276 633 F 1300 998 830

Customer Service New Zealand
T 0508 276 111 F 0508 998 830

©Registered Trademark

BSNmedical

THERAPIES. HAND IN HAND.

BSNmedical

LEUKOSAN® STRIP


Wound closure strips are an ideal option for topical applications, providing atraumatic, painless wound closure, low infection risk and good cosmetic results. Thanks to our knitted elastic fabric, **Leukosan® Strip** represents the latest state of development, enabling improved clinical performance.

RELIABLE ADHESION FOR LONGER WEAR TIME

The permeable fabric allows the skin to breathe and moisture to escape. The strong, yet skin friendly adhesive keeps the strips securely in place.

FRAY-PROTECTED FOR GREATER DURABILITY

Rounded edges prevent material from fraying or rolling up and assure durability of strips thus minimising aftercare.

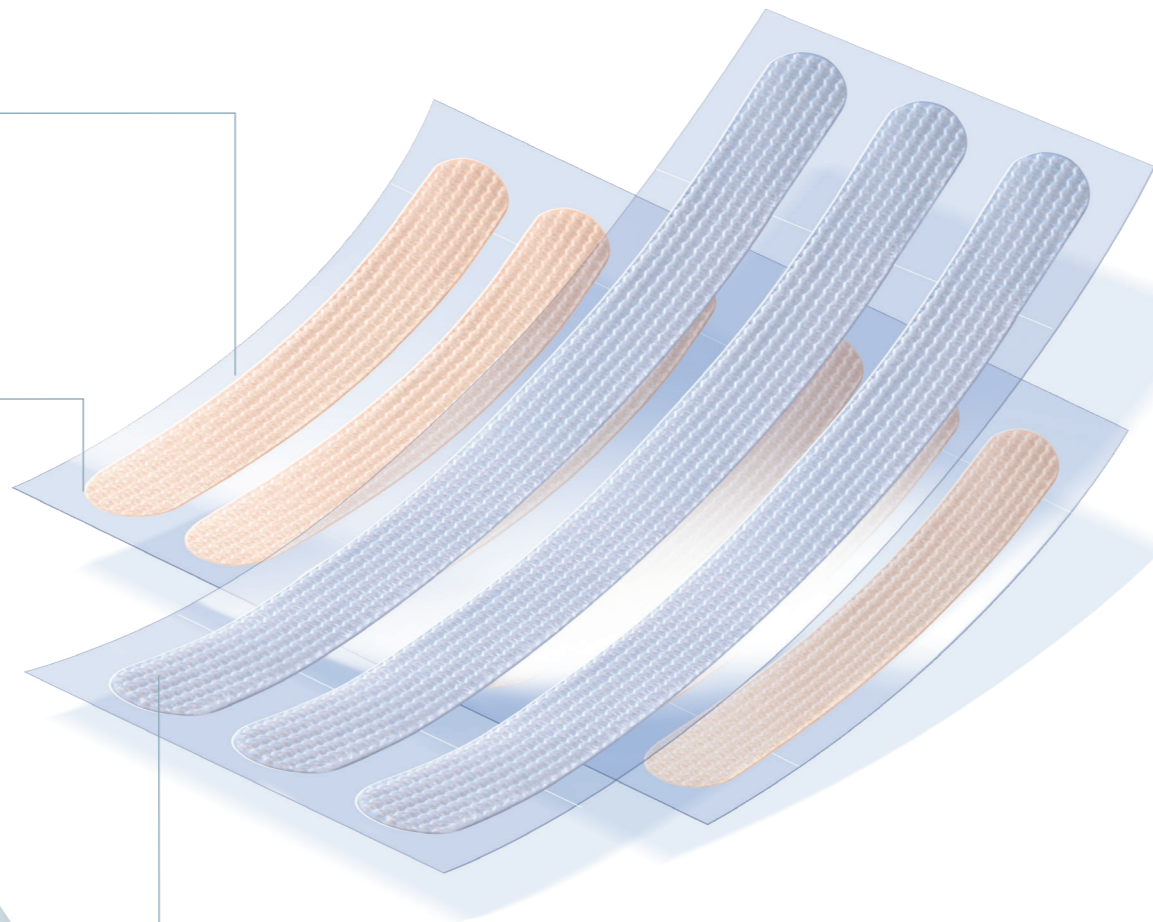


INTELLIGENT THREADS FOR ENHANCED COMFORT AND CONFORMABILITY

Curved polyamid threads enable strips to expand and retract. This feature allows strips to adapt to oedema-caused expansion of the wound tissue whilst avoiding constriction.



Also available in Tan colour

Leukosan® Strip Tan is less visible in exposed areas like the face and the hands.

Leukosan® Strip adapts to the changing contours of the wound as it heals (1). It expands with the oedema phase (2) and will return to its initial length when swelling subsides (3). As a result, the wound edges are kept securely aligned at all times, promoting more effective healing.

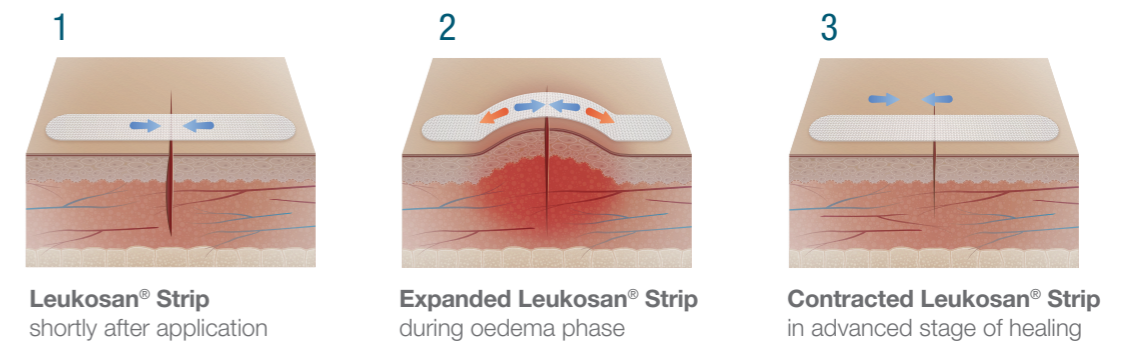

External forces caused by oedema
Retractive power of thread

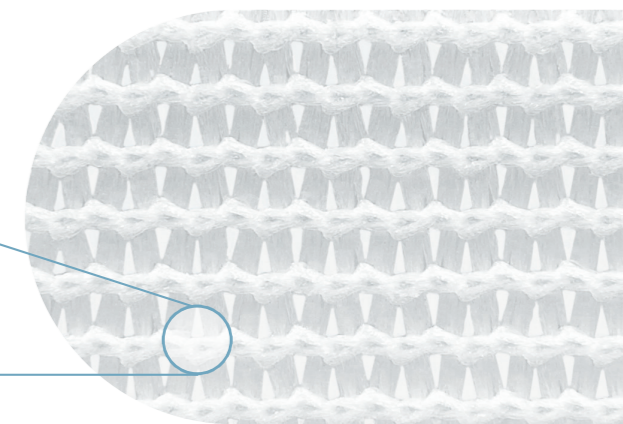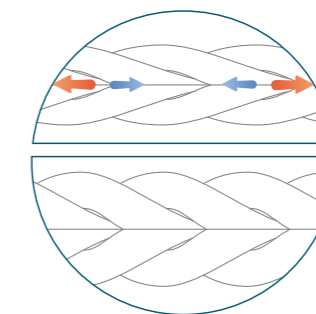

Leukosan® Strip is recommended for the following indications



Leukosan® Strip is ideal for closing facial lacerations or incisions without additional trauma and best cosmetic results.



Additional tension relief for wounds closed with intracutaneous sutures.



For use on thin, fragile or sensitive skin.



Early replacement of staples or sutures to improve cosmetic results.